



DGA Medical, LLC
A Medical Parts & Equipment Company
602-266-6389

A Veteran-Owned Small Business

Biodex Basic Ultrasound Table

Designed for all Ultrasound and OB/GYN procedures.

- Ergonomic design helps prevent musculoskeletal injuries
- Open frame for greater accessibility
- Contoured tabletop allows optimal probe positioning
- Features three primary positions: horizontal, chair and Trendelenburg
- Manual positioning for height, Fowler back and Trendelenburg motions
- Optional side rails fold down for unobstructed patient access

The Basic Ultrasound Table features an hourglass shaped tabletop that narrows to 19 inches, providing comfortable scanning probe positioning at and below the plane of the tabletop. The leg section unlocks and swings down for chair positioning, or for positioning with optional stirrups.

Table adjustments are straightforward on the Basic Ultrasound Table. To adjust table height, a simple hydraulic foot pump is used. Trendelenburg positioning is easily performed with a smooth-turning hand crank and Fowler positioning is achieved with preset stops at 0°, 15°, 35°, 55° and 80°.


Suitable Procedures:

- Abdomen
- Aorta
- Endovaginal*
- Hysterosonography*
- Testicular
- Musculoskeletal
- DVT/Lower Extremity Venous*
- Renal Artery
- Abdominal Doppler
- Transesophageal Echo
- Renal
- Pelvis
- OB
- Thyroid
- Breast
- Carotid
- Lower Extremity Arterial
- Renal Transplant
- SMA

*Requires optional stirrups.

Specifications

- Dimensions: 75" l x 24" w (191 x 61 cm); 30" w (76 cm) with side rails; 19" w (48 cm) at narrow section
-
-

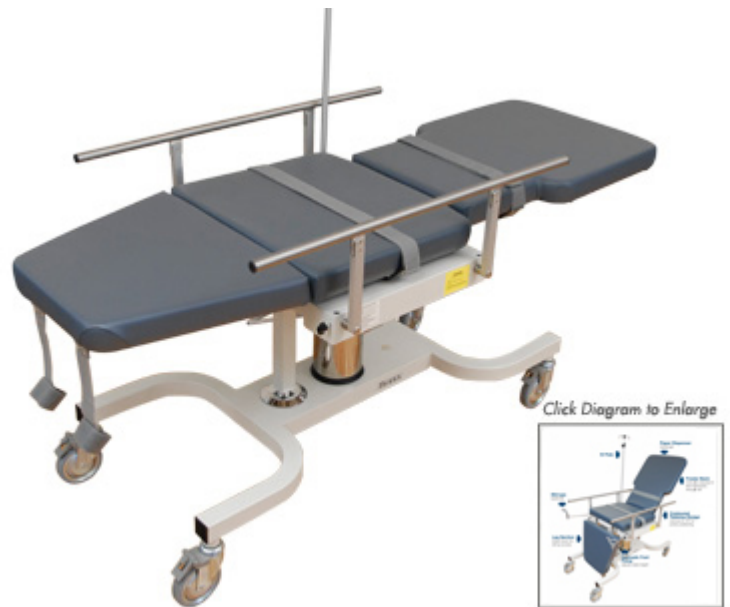

Biodex™ is a registered trademark of Biodex Medical Systems, Inc.
© Copyright Biodex Medical Systems, Inc. 2010

- Drop-Down Leg Section: 14.5" l (37 cm), drops down for stirrup access
- Height Adjustable: 29.5" to 37" (75 to 94 cm)
- Tabletop: Three primary sections; torso section features contoured design with Fowler back positioning, center section remains fixed and leg section drops down for stirrup access
- Mattress: 2" (5 cm) thick; foam rubber pad covered with durable Naugahyde®
- Upholstery Color: Graphite
- Motions:
 - Height Adjustable: 29.5" to 37" (75 to 94 cm), hydraulic foot pump
 - Trendelenburg: 0° to 15°, hand crank
 - Fowler Back: Manually adjusted in five preset increments, 0°, 15°, 35°, 55° and 80°
- Wheels: 5" (13 cm) locking casters
- Patient Restraints: Two body straps
- Finish: Powder coat
- Patient Capacity: 350 lb (159 kg); weight tested to four times the patient load rating.
- Weight: 335 lb (152 kg)
- Warranty: Two years



056-613 Table, Ultrasound, Basic with Convenience Package
Includes: Side Rails, Stirrups and IV Pole

056-609 Table, Ultrasound, Basic
Includes: Side Rails

056-605 Table, Ultrasound, Basic
*All options available separately.**

Optional:

056-608 Stirrups

056-607 IV Pole

058-610 Dispenser, Paper with Cutter Strap

Replacement:

056-612 Table Paper, Hygienic, 12 rolls
18" w x 225' l (46 cm x 69 m)

**Note: Options are installed by the customer.*

