



**DGA Medical, LLC**  
A Medical Parts & Equipment Company  
602-266-6389

## Biodex Echo Pro™ Echocardiography Table

### *Accommodates echocardiography procedures*

- Ergonomic design helps prevent musculoskeletal injuries – comfortable shoulder, arm and hand positions when scanning
- Allows seated or standing work positions
- Extra wide top with 500-lb patient capacity, suitable for bariatric patients
- Height range down to 23" accommodates wheelchair transfers and those with ambulatory difficulties
- Motorized height, Fowler back, Trendelenburg and auto level adjustment
- Swing-out cardiac scanning cutout doubles as armrest
- Adjustable sonographer cutout for maximum patient access and comfortable scanning positioning
- Vascular scanning arm board mounts to either side of table

There are no barriers between sonographer and patient on the Echo Pro table, allowing correct scanning posture to prevent musculoskeletal injury. The cardiac scanning cutout provides open access to the patient's left thorax area for an unobstructed apical approach. The motorized Fowler back is infinitely adjustable to 80° via hand or foot controller.

The Echo Pro Combination Table is so accommodating it reduces the total amount of time required to achieve a quality image. Patient comfort is assured even through lengthy procedures by cushioned arm supports, a thick mattress and a unique swing-out cardiac scanning cutout cushion that doubles as an armrest.


#### **Suitable Procedures:**

- Abdomen
- Aorta
- Thyroid
- Breast
- Carotid
- Venous Reflux\*
- Lower Extremity Arterial
- Upper Extremity Arterial\*\*
- Renal Transplant
- SMA
- Stress Echo
- Renal
- Pelvis
- Testicular
- Musculoskeletal
- DVT/Lower Extremity Venous\*
- Upper Extremity Venous
- Renal Artery
- Abdominal Doppler
- Echocardiography
- Transesophageal Echo

\*Requires optional foot support

\*\*Requires optional vascular arm board

#### **Specifications**

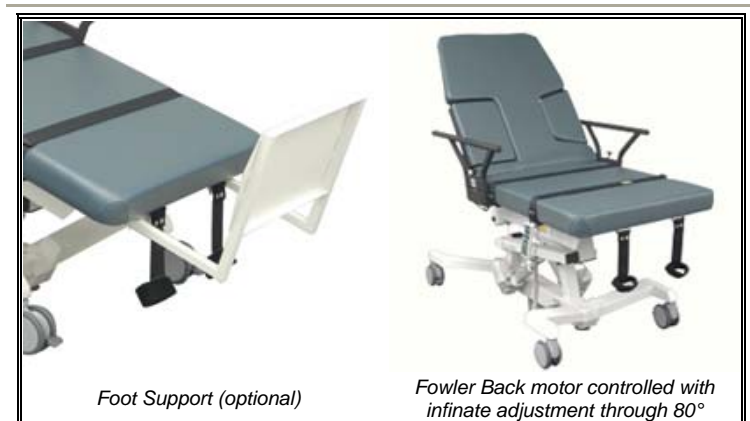

- Dimensions: 73" l x 30" w (185.4 x 76.2 cm); 35" w (88.9 cm) with optional side rails
- Height Adjustable: 23" to 39" (58.4 to 99 cm)
- Tabletop: Two primary sections: torso section features sonographer and patient cutouts with Fowler back positioning, and lower body section remains fixed
  - Cardiac Scanning Cutout: with pivoting cushion that doubles as an armrest
  - Adjustable Sonographer Cutout: drops down or folds up for maximum access to the patient
- Mattress: Torso section 2" (5 cm) thick, lower body section 3" (8 cm) thick; foam rubber pad covered with durable Naugahyde®
- Upholstery Color: Graphite
- Motions:
  - Height Adjustable: 23" to 39" (58.4 to 99 cm)
  - Trendelenburg: 0° to ±15°
  - Fowler Back: 0° to 80° infinitely adjustable
- Controls:
  - Hand Control: Activates height, Trendelenburg motions, Fowler back and auto level motions
  - Foot Control (optional): Activates height, Trendelenburg motions, and Fowler back
- Wheels: 5" (12.7 cm) Individual locking swivel casters; central floor locking system
- Patient Restraints: Two body straps
- Finish: Powder coat
- Patient Capacity: 500 lb (227 kg); weight tested to four times the patient load rating.
- Weight: 335 lb (152 kg)
- Power: 115 VAC or 230 VAC
- Warranty: Two years
- Certifications: ETL and cETL Listed to UL60601-1 and CAN/CSA C22.2 No. 601.1-M90 and EN60601-1 standards.




---

**058-700** Table, Echocardiography, Echo Pro™  
115 VAC  
Includes Hand Controller.  
*All options available separately\**

---

**058-705** Table, Echocardiography, Echo Pro™  
230 VAC  
Includes Hand Controller.  
*All options available separately\**

Optional:

---

**058-740** Foot Controller

---

**058-632** Side Rails, Folding

---

**058-735** Arm Board, Vascular Scanning

**058-** IV Pole  
**656**

**058-** Foot Support  
**640**

**058-** Dispenser, Paper  
**610** with Cutter Strap

Replacement:

---

**056-** Table Paper, Hygienic, 12 rolls  
**612** 18" w x 225' l (46 cm x 69 m)

---

*\*Note: Options are installed by the customer.*

US Patent 7,103,932 B1